Q-Temp™ / Q-Temp Clear™

Non-eugenol, self-cure resin based temporary cement

For professional dental use only

EN INDICATIONS

- Q-Temp is used for temporary cementation of crowns and bridges, inlays and onlays.

PROPERTIES

- Q-Temp is a self-curing cement.
- Q-Temp classification according to ISO 4049 is Type 2, Class 1.
- Ensures firm retention.
- Easy Removal of excess cement.
- Q-Temp has been on the market since 2008 with no changes in formula. No complaints of any health-related or chemical performance issues were received.

CONTRAINDICATIONS

- Q-Temp Temporary Cement is not intended as a temporary filling material nor is it to be used for pulp capping.
- Q-Temp Clear has a significantly low radiopacity due to low filler content.

CAUTION!

- Though rare, some individuals may be sensitive to this product. If an allergic reaction occurs, discontinue use, remove from patient and refer patient to a physician.
- Q-Temp contains materials which: irritating to eyes / irritating to respiratory system and skin. For complete symbols guide see below.
- Large amount of excess material may generate heat as part of the normal setting reaction.

COMPOSITION Q-TEMP

Essential ingredients – BASE	% wt.	Essential ingredients – CATALYST	% wt.
Urethane dimethacrylate	<40%	Urethane dimethacrylate	<50%
Polymerization activators	<5%	Benzoyl Peroxide	<5%
Fumed silica	<5%	Fumed Silica	<5%
Titanium oxide	<1%	Titanium oxide	<1%

FILLER CONTENT Q-TEMP

The filler system comprises modified silica and minerals:
10-20 weight-% or 2-10 volume-%.
The particle size of the fillers range: <6 µm.

COMPOSITION Q-TEMP CLEAR

Essential ingredients – BASE	% wt.	Essential ingredients – CATALYST	% wt.
Urethane dimethacrylate	<40%	Urethane dimethacrylate	<50%
Polymerization activators	<5%	Benzoyl Peroxide	<5%
Fumed silica	<5%	Fumed Silica	<5%

FILLER CONTENT Q-TEMP CLEAR


The filler system comprises modified silica and minerals:
Less than 5 weight-% or 2 volume-%.
The particle size of the fillers range: <6 µm.

INSTRUCTIONS FOR USE

Note:

Tips are for single patient use only to prevent cross-contamination between patients.

- Complete tooth preparation and fabrication of the temporary restoration before extruding the mixed cement.
- Dry and isolate the prepared tooth. Be careful not to desiccate or over-dry the preparation.

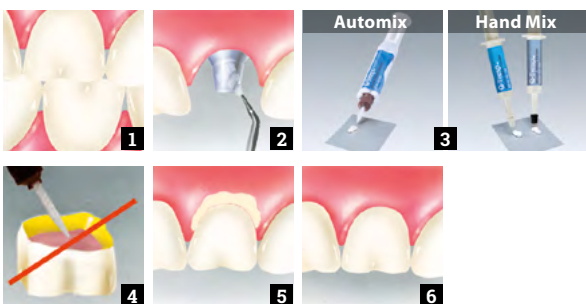- Automix (When using a dual-barrel syringe with automix syringe mix tip):** a) Bleed the syringe before installation of the automix syringe mix tip. Mount a new, automix syringe mix tip onto the dual-barrel syringe and bleed the syringe again. This ensures complete flow of Base and Catalyst components. b) Syringe a thin layer of Q-Temp on to the cleaned, inside surface of the temporary restoration.

- Hand Mix (When using a two syringes system of Base and Catalyst):** a) Extrude equal amount of Base and Catalyst onto the selected mixing surface. b) Using a clean spatula or mixing stick, mix the two components into a homogenous paste. c) Using the mixing stick or brush, apply a thin layer of the mixed Q-Temp on to the cleaned, inside surface of the restoration.

Note: Only a small bead of Q-Temp around the margin is necessary to securely hold your provisional. This is especially true with highly retentive temporaries. Using less cement also allows easy retrievals.

- Firmly set the temporary restoration in place until the initial set. The initial set is characterized by a distinct gel state. Excess cement is best removed in its gel state, after reaching the material working time, with a scaler or explorer.
- Tightly re-cap syringe after use.
- For the material time indications, see table below.

Time indications:

The working and setting times depend on the ambient and oral temperature.

The times shown are based on conditions relevant for practice.

	min:sec
Working time at ambient temperature	>01:30
Intra-oral setting time	<05:00

Mixing options:

Note: working and setting times complies to the mentioned above range without dependence on mixing options.

Mixing Options	Ratio
Automix - 1 Automix syringe	1:1 - Equal amount of Base and Catalyst
Hand Mix - 1 Automix syringe with Detachable Manual Adapter	1:1 - Equal amount of Base and Catalyst

STORAGE AND DISPOSAL

- Do not store above 24°C (75°F).
- Keep product in original packaging.
- For optimum freshness, keep refrigerated.
- Do not freeze.
- If refrigerated, allow the syringe to reach room temperature.
- To obtain SDS or IFU visit www.bjmlabs.com.
- Do not use after expiration date. See expiry date on the label.
- Shelf life is 2 years after the date of manufacturing.
- Empty product can be safely disposed after use in an intended regular waste container. The product does not require special disposal instructions or precautions and is harmless to the environment.

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